

C H E V A L I E R

Hydraulic STRETCHER



MEDEGEN
Medical Products
An Inteplast Group Company

FEATURES



Wave-Shape Siderail

Collapsible stainless siderail with ABS wave-shape handle is rust-proof and easy to grasp and clean.

Siderail Safety Release

Pull out the locking safety release before folding down siderails to ensure patient safety.



ACCESSORIES



Foot-End Shelf

- Available to set EKG monitor or become a writing board.
- Easy to store when not in use.



Hi-Lo Function

- Hydraulic system with sophisticated mechanism provides a steady, durable, safe and comfortable elevation function.
- The height of beds can adjust from 22.05° to 35.04° by stepping on the hi-low pedal located on both sides of stretcher.

Trend / Reverse Trend Function

The Trend/Reverse Trend function can adjust from -16 to 16 degrees with stability and safety at any height of the stretcher.



Chair Position Function

- Settable Chair Position.
- Various position choices offered to improve patient comfort.



1 I.V. Pole

Two-step adjustable stainless I.V. pole (4-section)

2 Foot Position Lever

Two-step foot lever enables foot section up to 15 degrees.



3 Name Card Holder (optional)

Located at the foot-end of the bed. Easy to place and replace patients' name cards.



4 Knee Position Crank

Manual crank is located at the foot-end and enables the knee section shifting from 0 to 40 degrees.



5 Base Cover

ABS Base Cover with storage bin.

6 Drainage Bag Hooks

Total 6 hooks available. 3 on each side of stretchers.



7 Foot Pedal

The foot pedal controls the height and tilt positions of stretcher by hydraulic mechanism.



8 Restraint Strap Holder

Total 10 holders available, 5 on each side of the stretcher.

9 Fifth Wheel (optional)

Fifth Wheel can improve the maneuverability of the stretcher.



Hydraulic STRETCHER

CHEVALIER

Model: CTH

10 Caster

8" central brake caster with dust-proof cover.

11 Central Braking and Steering Pedal

Enables stretchers to move, brake and steer.



14 Back Position Handle

The backrest is easily positioned from 0 to 90 degrees using a dual gas-spring set below the backrest platform.



15 Driving Handle

Driving handle can be located at either head, foot or both ends for smooth movement and can be stored while not in use.



16 ABS Corner Bumper

4 rounded ABS corner bumpers are set at each corner of the stretcher.

17 Oxygen Tank Holder (optional)

0.5 m³ Oxygen Tank Holder, located at the head section, can rotate in/out, according to medical requirement.



12 Side Rail


ABS Collapsible Siderail

13 X-ray Platform


The X-ray platform is available for medical treatment in two choices:

1. Movable full-length X-ray platform.
2. Compact backrest X-ray platform.




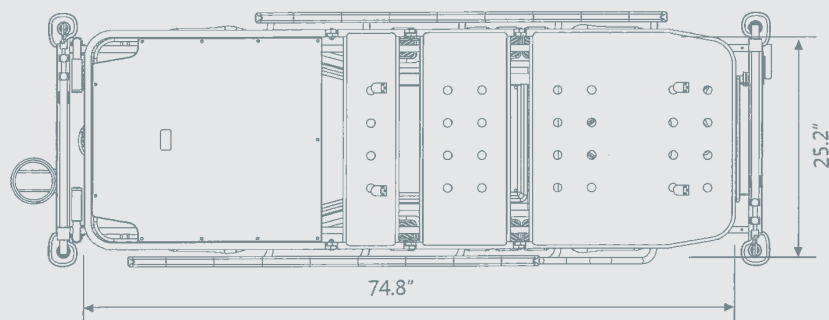

Top View

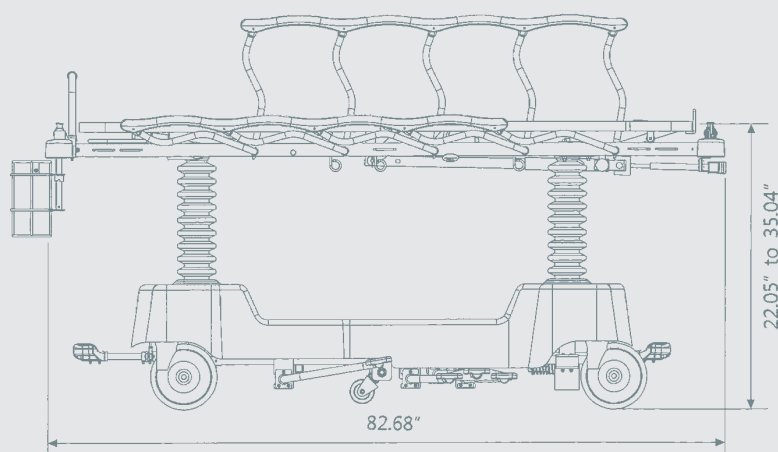

Front View

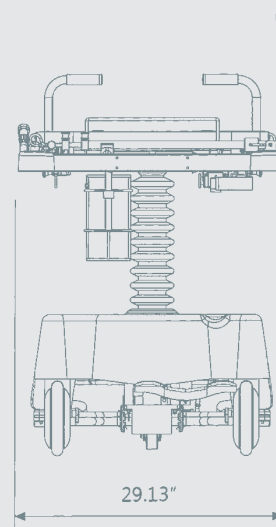

Right View

SPECIFICATIONS

Dimension

Overall Size	L 82.68" x W 29.13"
Sleeping Surface	L 74.8" x W 25.2"
Siderail	L 47.24" x W 14.57"

Positioning

Hi-Low Range	22.05" – 35.04"
Backrest	0 – 90°
Knee Gatch	0 – 40°
Trend / Reverse Trend	0 ± 16°
Max. Weight Capacity	230 kgs (500 lbs)

Standard Accessories

Driving Handle	1ST
Fixed I.V. Pole	2PC

Optional Accessories

Mattress
Mattress Cover
Fifth Wheel
0.5 m ³ Oxygen Tank Holder
Name Card Holder
Foot-End Shelf

LIMITED 2 YEAR WARRANTY



209 Medegen Drive, Gallaway, TN 38036
Toll Free 800-511-6298
Customer Service Fax 631-714-6311
E-mail cservice@medegenmed.com



INTEPLAST GROUP
Healthcare

© 2022 Medegen Medical Products LIT# 0076B