

Biological COLLECTION

Starswab II™ collection is a complete bacteriology and transport system. The color coded cap has a protective over skirt which helps ensure the containment of biohazardous material.

- All Starswab II™ packaging is labeled with product number, lot # and expiration date.
- The vials are made of polypropylene with a thick tube wall to protect medium from dehydration and the spiral media chamber keeps sample immersed and protected during transport.
- The round bottom tube makes for easy racking and the seamless tube molding ensures tube strength and integrity.
- The rayon tip is inert, non-toxic and permits good sample retrieval and absorption.
- Amies Clear provides excellent survival for a wide range of bacteria and is compatible with gram stain.
- Amies Clear provides excellent environment for the survival of anaerobes through the complete immersion of the swabs into the gel and by the presence of the reducing agent, sodium thioglycollate
- Modified Liquid Stuart's and Liquid Amies provide excellent recovery in routine bacteriology and are compatible with gram stain and many rapid screening methods.

Starswab II™ packaging promotes optimal medium performance. A premium thick bag creates a barrier environment to prevent premature medium dehydration. The metalized bag further protects the medium from harmful UV rays. The bag also features a ziplock closure which simplifies resealing to better maintain the remaining swabs at peak performance.

| Code # | Medium | Swab | Cap Color | Sterile | Case Pack | Packaging |
|----------------------|-----------------------------|--|---------------|----------|------------|---------------------------|
| P02-SP130X | Amies Clear | Single, Plastics Shaft, Rayon | Magenta | Y | 500 | 10 Bags of 50 each |
| P02-SP130D | Amies Clear | Double, Plastic Shaft, Rayon | White | Y | 500 | 10 Bags of 50 each |
| *P02-SP130L-D | Amies Clear w/sponge | Double, Plastic Shaft, Rayon | White | Y | 500 | 10 Bags of 50 each |
| P02-SP132 | Liquid Stuart's | Single, Plastic Shaft, Rayon | Teal | Y | 500 | 10 Bags of 50 each |
| *P02-SP132D | Liquid Stuart's | Double, Plastic Shaft, Rayon | Yellow | Y | 500 | 10 Bags of 50 each |
| *P02-SP132-W | Liquid Stuart's | Single, Straight Wire Shaft, Fine Tip | Orange | Y | 500 | 10 Bags of 50 each |

*Item is new Medegen Medical Product offering

MEDEGEN
Medical Products
An Inteplast Group Company

209 Medegen Drive, Gallaway, TN 38036
Toll Free 800-511-6298
Customer Service Fax 631-714-6311
E-mail cservice@medegenmed.com




INTEPLAST GROUP
Inteplast Healthcare



P02-SP130L-D



P02-SP132-W

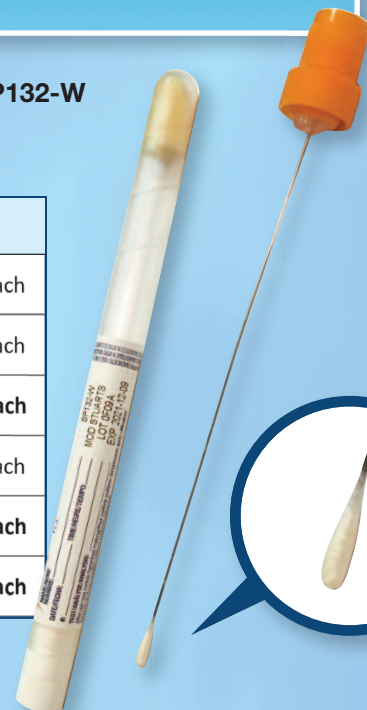

Biological COLLECTION

Medegen Medical Products has expanded their Biological collection offering.

- Multitrans media with glass beads is packaged in a durable Polypropylene vial with a skirted swab retaining cap.
- Swabs are scored so it is easily snapped off into the vial; after capping and reopening, the swab is neatly inserted into the cap for easy, contamination free handling.
- Flocked swabs ensure superior specimen collection and release which provides a more representative specimen sample for testing.
- Flocked swab technology optimizes patient comfort during specimen collection and provides a more representative sample for testing.
- The media vial and sterile wrapped swab are packaged in a sterile plastic peel pack.



SW-030FS




P02-S160-100



P02-S160-NASO



P02-S160-0

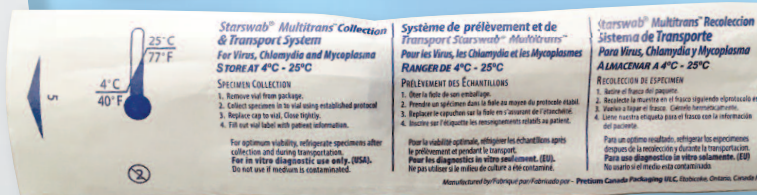

SW-080FS



| Code # | Media | Swab | Cap Color | Sterile | Case Pack | Packaging |
|---------------|------------|------------------------------------|-----------|---------|-----------|-------------------|
| P02-S160-100 | Multitrans | Single, Plastic Shaft, Polyester | Blue | Y | 100 | 2 Bags of 50 each |
| P02-S160-NASO | Multitrans | Single, Plastic Shaft, Fine Tipped | Blue | Y | 100 | 2 Bags of 50 each |
| P02-S160-0 | Multitrans | None | Blue | Y | 100 | 2 Bags of 50 each |

P02-S160-0 is media only. We also offer sterile wrapped flocked swabs.

| Code # | Description | Sterile | Case Pack | Packaging |
|----------|--------------------|---------|-----------|----------------|
| SW-030FS | Flocked Oral Swab | Y | 1000 | 5 boxes of 100 |
| SW-080FS | Flocked Nasal Swab | Y | 1000 | 5 Boxes of 100 |

Swabs are enlarged to show detail.



209 Medegen Drive, Gallaway, TN 38036
Toll Free 800-511-6298
Customer Service Fax 631-714-6311
E-mail cservice@medegenmed.com



INTEPLAST GROUP
Inteplast Healthcare

