


C H E V A L I E R

Hydraulic STRETCHER



 **MEDEGEN**
Medical Products
An Inteplast Group Company

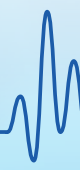

Trend / Reverse Trend Function

The Trend/Reverse Trend function can adjust from -16 to 16 degrees with stability and safety at any height of the stretcher.

Hi-Lo Function

- Hydraulic system with sophisticated mechanism provides a steady, durable, safe and comfortable elevation function.
- The height of beds can adjust from 22.05° to 35.04° by stepping on the hi-low pedal located on both sides of stretcher.



Chair Position Function

- Settable Chair Position*.
- Various position choices offered to improve patient comfort.

*Only available on stretcher with ½ body radiolucent board

Wave-Shape Siderail

Collapsible stainless siderail with ABS wave-shape handle is rust-proof and easy to grasp and clean.



Siderail Safety Release

Pull out the locking safety release before folding down siderails to ensure patient safety.

1 I.V. Pole

Two-step adjustable stainless I.V. pole (4-section)

2 Foot Position Lever

Two-step foot lever enables foot section up to 15 degrees.



3 Built-in Scale

Comes standard on stretcher CSNN.



4 Knee Position Crank

Manual crank is located at the foot-end and enables the knee section shifting from 0 to 40 degrees.



5 Base Cover

ABS Base Cover with storage bin.

6 Drainage Bag Hooks

Total 6 hooks available. 3 on each side of stretchers.

7 Foot Pedal

The foot pedal controls the height and tilt positions of stretcher by hydraulic mechanism.



8

8 Restraint Strap Holder

Total 10 holders available, 5 on each side of the stretcher.

9 Fifth Wheel

Fifth Wheel can improve the maneuverability of the stretcher.



10 Caster

8" central brake caster with dust-proof cover.

12

11 Central Braking and Steering Pedal

Enables stretchers to move, brake and steer.



12 Side Rail

ABS Collapsible Siderail

13 X-ray Platform

The X-ray platform is available for medical treatment in two choices:

1. Movable full-length X-ray platform.
2. Compact backrest X-ray platform.

14 Back Position Handle

The backrest is easily positioned from 0 to 90 degrees using a dual gas-spring set below the backrest platform.



15 Driving Handle

Driving handle can be located at either head, foot or both ends for smooth movement and can be stored while not in use.



16 ABS Corner Bumper

4 rounded ABS corner bumpers are set at each corner of the stretcher.

17 Oxygen Tank Holder (optional)

0.5 m³ Oxygen Tank Holder, located at the head section, can rotate in/out, according to medical requirement.



Med Stretch Plus

REPLACEMENT MATTRESSES



GMT4TVNE

Premium Mattress 74.8" X 25.2" X 5"

- Material PU
- ANTI-BACTERIAL
- 2 layers of foam including soft comfort memory foam for patient well-being
- Easy to Clean
- Latex Free
- For Medegen Stretchers HYDSTRETCHER and CSNNSTRETCHER



***Premium replacement mattresses
for industry standard stretchers***

GMT4TW5A


Premium Mattress 79" x 26" x 5"

GMT4TX5A


Premium Mattress 79" x 29" x 5"

- 2 layers of foam including soft comfort memory foam for patient well-being
- Material PU
- ANTI-BACTERIAL
- Easy to Clean
- Latex Free
- Mattress break point reducing pressure points and skin break down.




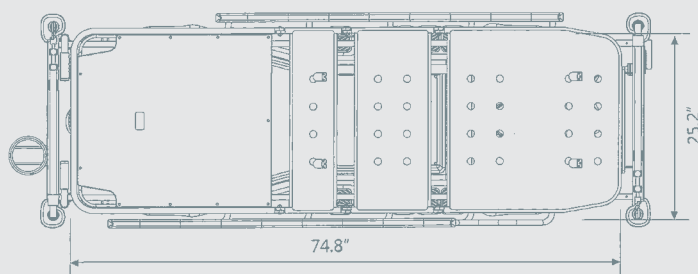

Top View

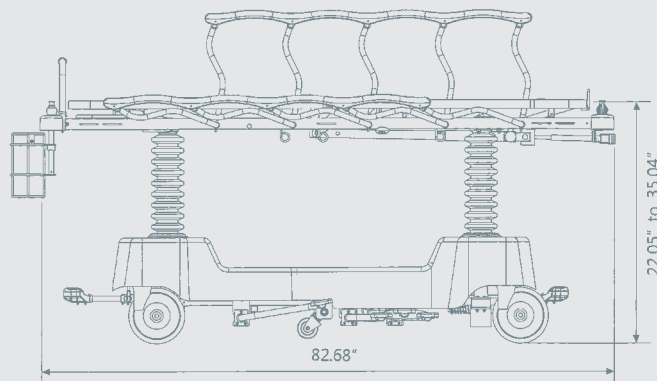

Front View

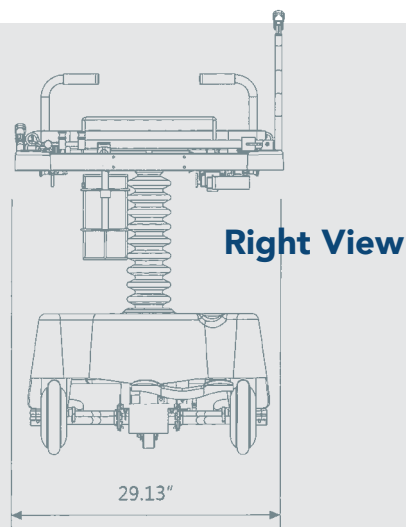

Right View

The ER stretchers we purchased from Medegen are incredible. They function very smoothly and the mattress is comfortable. The flat surfaces make them easy to clean for a quick turnaround. Our rep came in and assembled the stretchers for us and gave a demonstration which is not very common out here in rural Texas. We appreciate the awesome service we received and look forward to continuing our relationship with Medegen.

Valerie Sanchez, Schleicher County Medical Center
102 N US Highway 277, Eldorado, TX 76936
325-853-2507 ext. 3929
Pharmacy / Central Supply-Purchasing / Housekeeping

SPECIFICATIONS

Dimension

Overall Size	L 82.68" x W 29.13"
Sleeping Surface	L 74.8" x W 25.2"
Siderail	L 47.24" x W 14.57"

Positioning

Hi-Low Range	22.05" – 35.04"
Backrest	0 – 90°
Knee Gatch	0 – 40°
Trend / Reverse Trend	0 ± 16°
Max. Weight Capacity	230 kgs (500 lbs)

Standard Accessories

Driving Handle	1ST
Fixed I.V. Pole	2PC

Optional Accessories

Mattress
Mattress Cover
X-Ray Cassette Holder
Paper Roll Dispenser
Foot-End Shelf
Stretcher Tray

LIMITED 2 YEAR WARRANTY



To view our products online, please visit us at www.medegenmed.com/carts or directly with your device by scanning the QR code to the right.



209 Medegen Drive, Gallaway, TN 38036
Toll Free 800-511-6298 • **Customer Service Fax** 631-714-6311
E-mail cservice@medegenmed.com

©2026 Medegen Medical Products LIT# 0076D